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Corporate Philosophy

This is what our name stands for

Modern life presents us with extraordinary challenges. Maintaining a healthy lifestyle allows us to face those challenges head on. In order to achieve optimal health we need to sleep well, get adequate exercise, and eat a healthy diet. Many people are afflicted with sleep disturbances and feel exhausted in the morning. Lack of sleep can adversely affect us both physically and emotionally.

It is not in vain that people say:

„Bed, or rest, is the best medicine.“

Eating the wrong foods can result in frequent illnesses or reduced energy. Lack of adequate exercise can lead to obesity, reduced stamina and disease. Sleep, exercise and diet, the three essential components of a healthy lifestyle, are the speciality of our company

Company founder Mr. Paul Neumann and his wife Mrs. Irene Neumann operate their business with their children „Star Naturwollprodukte und mehr...“ e.K.“ headquarters in Schwäbisch Gmünd-Bargau.

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Since 2000 we as a company appreciate the trust that is placed in us and our products.

We currently employ several hundred workers in our subsidiaries in Germany.

We also have subsidiaries in Austria, Italy, Belgium, Poland, the Czech Republic, Greece, Ukraine, Russia, Israel and other international locations.

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Royal-Deluxe Massage Bed

(Height: 29 cm)

with oscillation and Shiatsu massage system

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Shiatsu, which originated in Japan, promotes relaxation and a sense of well-being. It enhances both recovery and regeneration of organisms. Shiatsu makes use of a variety of pressure techniques applied to fingers, hands, or elbows. In Shiatsu, the knowledge of traditional techniques is combined with modern methods. We have put together all of these pressure techniques in one complete mechanical unit, so that anyone can utilize this relaxing practice. In our star-medico® Royal-Deluxe Massage Bed you have the options and advantages of oscillation therapy combined with an oasis for your personal needs.

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Ergoflex Comfort 4-pieces

(Height: 21 cm or 29 cm)

with oscillation (acknowledged Therapy Method)

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A new age of healthy sleep has begun. Our company recognizes the relationship between optimal sleep and well-being by combining hygiene, comfort and therapy. Your body will thank you when you sleep well, and you will be rewarded with clarity and energy when you awaken. The integrated special oscillation for common health complaints can be easily operated utilizing pre-installed programs. You have a choice whether you want to stimulate only certain areas of the body or the whole body with pleasant, relaxing vibrations.

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Tencel®

Botanical fiber from eucalyptus wood

Tencel® has a unique way of controlling humidity because of the fine nanofibril in the internal structure of the threads. These threads act as a kind of microscopic canal system which transports the body's humidity away from the skin making you feel comfortable and relaxed

In addition, this unique system **Silver Vital Plus®** ensures that you will be protected from mildew, bacteria and mites while also benefiting from Tencel's odor-free design without the addition of chemicals.

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Perfect climate qualities

STEAM ADMISSION WITH EXTREME AIR HUMIDITY

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SilverVital Plus®

“Silver is worth gold not only during sleep”

Silver Vital Plus® is our own brand, which is a very high quality natural fiber with 3% Silver thread share. Until the beginning of this century silver was recognized for its medicinal properties. Today, its healing properties are being rediscovered. Within four to six minutes, colloidal silver eliminates a variety of fungal infections, viruses, bacteria, streptococcus, staphylococcus, and other pathogens such as house dust mites. It also serves as an odor neutralizer. According to current scientific knowledge there is no known bacteria which cannot be eliminated by colloidal silver. Even in high concentrations it has no side-effects.

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star-medico[®] Ergoflex Anti-Allergenic Cover

The clever three-layer model

1. The first outside layer is made of an elastic and pre-washed Tencel[®] fabric with Silver Vital Plus[®], that does not shrink when washed.

2. The second layer is made of a uniquely developed hollow, thermal, and breathable fiber which is heat-insulated. The fiber stores body heat, and the stored heat is maintained, keeping you warm. At the same time its distinctive air circulation system inhibits unpleasant perspiration.

3. The third internal layer is composed of a high-quality functional fabric. It is extremely point-elastic, breathable, and antibacterial which due to cost, is not generally used in convention mattress covers.

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star-medico® Ergoflex mattress

with pressure relieving slat technology

The exclusive star-medico® Ergoflex is a high point-elastic orthopedic mattress, which offers full body support. With the star-medico® Ergoflex, heavier parts of the body such as the shoulders and pelvis can easily sink into the core, while the lighter body parts like the lumbar spine and the knee joint can be supported individually. With this level of support and the level of relaxation that it promotes, pressure on the lumbar spine is much reduced. As a consequence, circulation improves and intra-vertebral discs are rejuvenated.

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star-medico® Pillow & Duvet

with the ThermoActive Tencel®-textiles
and Silver Vital Plus®

ThermoActive-pillow:

- Height can be individually adjusted through the stitching and filling

OrthoRelaxe:

- Pillow visco foam

ThermoActive-pillow and duvet :

- Washable at 60°C
- Antiallergic and antistatic
- Material contains 3% silver thread
- Tencel® and Silver Vital Plus® climate regulation
- Adapts perfectly to bed shape
- Cover can be washed separately
- Excellent durability
- Size: 135/155/180/200 x 200/220/240 cm
- 5 year warranty
- Quality: Made by Star Naturwollprodukte und mehr...*e.K.

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star-medico® Pad Flex Ortholine frame

Ergonomic reclining comfort

An orthopedic mattress should have a well-functioning bed insert to ensure optimal performance. The flexible star-medico® Pad Flex Ortholine frame serves as a perfect base for the orthopedic mattress and should be necessarily used for optimal storage. In addition to our Sleep-Therapy-System (STS), the star-medico® Pad Flex Ortholine frame, is equipped with a manual or motorized adjustment.

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2. Ergoflex Duvet Cover

3. "Bioactive"
silver reference

4. Upper bed
"ThermoActive Silver
Vital Plus"

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Oscillation (accredited therapy method)

Advantages of oscillation

Our orthopedic sleeping pad with oscillation (STS) can help with several complaints:

- It provides good relaxation, even during a short sleep,
- it helps to initiate sleep, and ensure uninterrupted sleep,
- it encourages the positive flow of lymphatic fluid,
- it reduces water retention in the legs,
- it firms-up tissue, and is helpful against cellulite,
- it stimulates circulation,
- it relieves joint pain,
- it relieves pain in the spine, and sciatica,
- it supports nutritional flow to the inter-vertebral discs,
- it promotes a sense of well-being,
- it increases concentration through better relaxation.

Due to the astounding success of the therapy, this method has found its way into the fields of health care and prevention.

- In fasting and as weight reduction
- In cleansing therapies
- In common regeneration
- As relaxation therapy
- In deacidification

Oscillation therapy is aimed at support the self-healing of the entire body, and has no negative side-effects.

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Infrared-Deep Heat

Regeneration and Anti-Stress-Action

Infrared heat stimulates blood circulation and expands blood vessels, even in non-radiated areas of the body. Blood pressure decreases, an cognition performance is improved. With blood flow, metabolic waste and acids are removed. The muscles become more relaxed, and tension pains are released.

The radiant energy of infrared heat allows your body temperature to increase without undue stress on circulation. Due to penetration of the subcutaneous tissue the necessary temperature can be reached. With an increase in body temperature, body fat is mobilized and toxins are expelled through the skin.

- Feel what helps -
Infrared abdominal massage belt
& Infrared heating pads

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star-medico® professional Massage Bed

with oscillation and
Infrared-Deep Heat

The regular use of an infrared heat system with oscillation therapy is an excellent way to detoxify and purify the body. Excessive acids, as fat-soluble environmental poisons, and heavy metals are secreted through the skin. By the use of Infrared deep heat the main constituents of perspiration are water, fat and cholesterol, including substantial amounts of heavy metals. Thus it was discovered that by using the infrared deep heat the amount of perspiration is doubled. As in addition the waterless part is tripled, the detoxifying potential can be up to six times more effective.

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Take up the fight against Cellulite

Many people suffer from it

What are the effects of oscillations on cellulite?

Oscillation therapy has a positive impact on all of the cellular and molecular processes of cellulite. It improves the flow of blood and lymph on all skin layers, and reduces infection and edema. The technical term is interrupted and the estrogen sensitivity of the skin cells decreases. It is proved that oscillation therapy is an effective treatment to reduce cellulite: The perimeter of the thigh and gluteus shrinks, the elasticity of the skin increases, and fibrosis is reduced.

Skin cross section
with cellulite

Skin cross section
without cellulite

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2005 Europa Fachschule Fresenius
Prof. Dr. Felder

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Implementation and studies

Professor Dr. habil. Vladimir Nazarov regarded as the founding father of the oscillation therapy based on BMS (biomechanical stimulation).

Throughout Germany many doctors' offices, rehabilitation, and health centers are already treating patients using the oscillation therapy method. Numerous studies prove that after a few weeks of use, a large number of patients notice considerable relief and are free of pain.

Studies* of supporting the effectiveness oscillation therapy:

2007 Clinic at the house lake Feldberg
Chief Physician Dr. Bucka MD,
Neurology, Chief Physician
Dr. Reich MD, orthopaedist

2005 Institute for Sport and Sports Science in the University of Karlsruhe (TH) Prof. Dr. Woll

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Contraindications

Oscillation therapy with Ergoflex Comfort is contraindicated to:

- persons with a permanent pacemaker, with acute thrombosis;
- persons after organ transplantation, in the presence of sepsis (blood poisoning);
- during pregnancy (recommendation: highly recommended after pregnancy);
- in severe bleeding (recommendation: during menstruation, the use of therapy is highly desirable);

Oscillation therapy program function - (OP) button 16

Press the (OP) button to activate treatment programs that can help with the following diseases.

OP1 Press the button (OP) – 1 time - for problems in the cervical and thoracic region of the spine, to relieve muscle tension and muscle tightness, in case of headaches and migraines.

OP2 Press the button (OP) – 2 times - in case of edema, fluid accumulation in the body, in order to reduce weight, to activate the lymphatic system, to treat hydroliipodystrophy (cellulite), in case of circulatory disorders.

OP3 Press the button (OP) - 3 times - for pain in the sciatic nerve, for problems in the lumbar spine, in case of lumbago, or diseases of the spinal column and spinal discs.

OP4 Press the button (OP) - 4 times - in case of arthritis, osteoarthritis, rheumatic diseases, joint pain, or osteoporosis; for metabolism stimulation.

OP5 Press the button (OP) - 5 times - during heavy physical exertion. You can use it for regenerating muscle pain, and for removing metabolic products from the body.

Functions of the health massage program. Button 12

P1 Press the button (P) - 1 time - for deep relaxation after physical stress;

P2 Press the button (P) - 2 times - to remove metabolic products, toxins, and toxic substances from the body;

P3 Press the button (P) - 3 times - to relax the muscles, in case of sleep disturbance, to restore strength;

P4 Press the button (P) - 4 times - for the energy and relaxation program. It has a supportive effect on performance, supports the overall process of regeneration and well-being, and stimulates the work of lymphatics very well.

P5 Press the button (P) - 5 times - for depressive symptoms, general malaise, or painful tension.


Heat source function. Button 13

Press the heat source button (W) to the ON position - or OFF (enable or disable). It provides additional relief from back pain and muscle tension. You can regulate the three-stage temperature mode by repeated pressing. Alternatively, switching using a separately sold belt is possible. Attention! Here only a belt or only a mattress is used. **Attention!** You can use only a belt or mattress.

Setup example

INSTRUCTION FOR EXPLOITATION

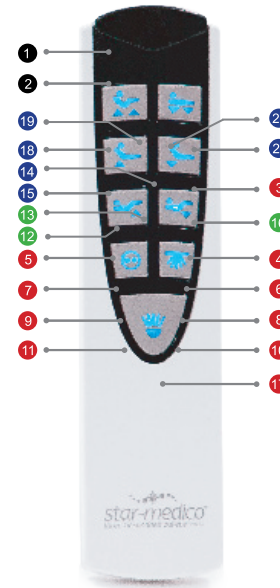- 1 Inspection LED indicator
- 2 LC-Display (LCD display)
- 3 "Massage" incl. / off
- 4 5 Massage module "upper back" (5) enable / increase intensity / intensity decrease / turn off
- 6 7 Massage module "lower back" (6) turn on/intensity (7) intensity decrease / turn off
- 8 9 Massage module for pelvis and thighs (8) turn on / increase intensity (9) intensity decrease / turn off
- 10 11 Massage module "lower legs. Feet" (10) turn on / increase intensity (11) intensity decrease / turn off
- 12 Wellness and massage programs 5 different programs
- 13 Thermal source. Three-stage temperature mode adjustable by repeated pressing
- 14 15 Wellness massage speed level (14) speed increase (15) speed reduce
- 16 Oscillating programs - 5 treatment programs
- 17 Light source on/off Timer "-"
- 19 sets the massage duration to 15 minutes (basic setting)
- 20 Timer "+" sets the duration of the massage to 30 minutes

Buttons for programmable functions:

- 18 Button - MENU (Menu) start programming mode
- 19 Button "-" changes the indicator down
- 20 Button "+" changes the indicator up
- 21 The PROG button advances to the next set-point or exits the programming mode

Made in Germany 

Riga: +37126371637
www.starmedico.lv

WE COMBINE HEALTHY SLEEP WITH THERAPEUTIC MASSAGE*****

Products are classified according to the German industrial standard for the reliability of electrical appliances for household use EN 60335-1. The use of prescribed standards guarantees compliance with the EC Low Voltage Directive and the electromagnetic compatibility of the devices. CE marking confirms the compliance with these standards. In addition, our products are subject to the regular controls. So far, no negative drug interactions have been observed.

Programmable Features

The following functions can be programmed:

- Time
- Alarm (beep)

For programming, only the buttons displayed at the top are used.


Time setting

- Press the **MENU** button. On the LC-Display (LCD) - the symbol "clock" flashes at the bottom left

- Press the **PROG** (program) button. On the LC-Display (LCD), the symbol "hours" and the numbers of the hour display flash.
- by pressing the + and – buttons on the time display set the desired time.

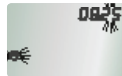- Press the **PROG** (program) button. On the LC-Display (LCD), the symbol "hours" and the numbers of the minute display flash
- by pressing the + and – buttons on the minute board set the desired time.

- Press the **PROG** (program) button
- The data is fixed, and the set time is shown on the display, the "clock" symbol disappears.

Disabling the alarm

- Press the **MENU** button and then the + button 1 time.
- The symbol "message" flashes in the middle left of the display

- Press the **PROG** (program) button
- the display shows **ON** (beep is on). Turn off the alarm with the **OFF** button (the sound signal turns off).

- Press the **PROG** (program) button. The data is fixed, and the display again shows the current time, the "message" symbol disappears.

Find us: „starmedico“

Program the sound signal (timer)

- Press the **MENU** button and then the + button 1x. The "message" symbol flashes on the display.

- Press the **PROG** (program) button
- The message symbol flashes on the display, the message **OFF** (beep off) and digits of the message time sensor.

- Turn the sound function on with the + button to **ON** (beep on).

- Press the **PROG** (program) button. The symbol flashes on the display "message" and message programs to choose from.


- Use the + and - buttons to set your preferred message program (User 1, Wave1, Wave 2, Wave 3, Pulse1, Pulse 2, or Pulse 3).


- Press the **PROG** (program) button.
- The "message" symbol and the hour display flash on the display.

- by pressing the + and – buttons on the time display, set the desired time.

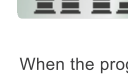- Press the **PROG** (program) button. The symbol flashes on the display. "message" and numbers of the minute board.
- by pressing the + and – buttons on the minute board, set the desired time.

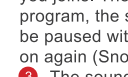- Press the **PROG** (program) button. Data is captured, and the display shows the current time, the time of the sound signal and the symbol "message".

When the programmed time of a sound signal approaches, the message program set by you joins. The program lasts 8 minutes, then it automatically turns off. During the message program, the symbol appears on the display "message". The message process can be paused with the button **17**. After about 5 minutes, the message program is switched on again (Snooze-Function). To turn off the sound function completely, press the button **3**. The sound function is activated the next day at the set time.

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Quality Control

“Confidence is good – Control is better!”

We collaborate with a variety of experts including professionals, doctors, therapists, sleep researchers and our suppliers. With their collaboration we continually develop and optimize our products.

Our products are produced according to the standards of TÜV SÜD, CE 1014, Textiles Vertrauen, TÜV Rheinland etc. (see below)

Under normal use, application, and by following the care instructions, you will receive a warranty lasting up to 5 years.

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Star-medico representative in the Baltics

Jurij Zubov

Doctor of natural medicine and massage.

Medical nutritionist. BMS (biomechanical muscle stimulation) therapist.

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